

-25°C Biomedical Freezer

DW-25L262



Hydrocarbon (HC) Energy-saving

DW-25L262 adopts an efficient hydrocarbon (HC) compressor and optimized refrigeration system, with power consumption as low as 0.95KW·h/24h.



Environmentally Friendly

The HC system is a completely green refrigeration solution. With LBA foamed polyurethane insulation in line with European ROHS directive for environmentally friendly materials.



Low Noise

Optimized hydrocarbon system and structural noise reduction design, DW-25L262 noise 34 dB, super quiet and low noise.



Security Lock

Mechanical lock plus lock catch, the mechanical door lock can be opened with one key, and the door lock can be padlocked. Secure and safe, suitable for multiple users.

-25°C Biomedical Freezer DW-25L262

Product Features



Microprocessor

- Microprocessor electronic temperature controller, LCD digital temperature display, display precision is 0.1°C, interior temperature set point -10~-25°C.
- Multiple alarms including high temperature alarm, low temperature alarm, sensor error alarm, power failure alarm, door ajar alarm and high ambient alarm.
- Alarm mode: sound and light alarm.
- It can be connected to remote alarm.



Hydrocarbon System

Increasing cooling rate, reducing energy consumption and lowering noise





Porthole

There is a standard access port (25mm inside diameter) on the left side of the unit, which is convenient for users to test independently.



Multiple Interface Options

- Standard configuration: Remote alarm interface.
- Optional configuration: USB or RS485.

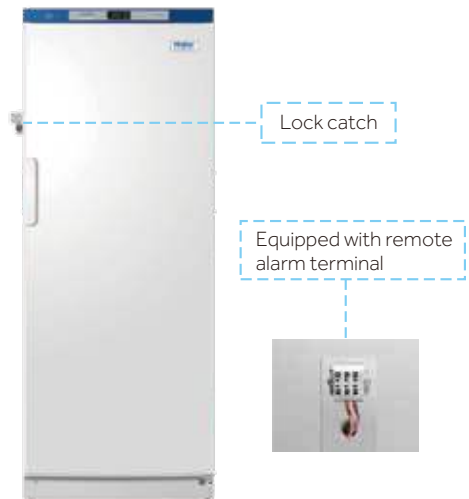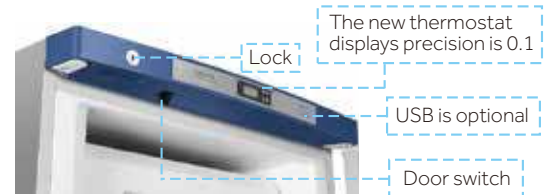

Key									
CODE	W0	W1	W2	W3	D0	D1	D2	H0	H1
MODEL									
DW-25L262	702mm	657mm	475mm	450mm	730mm	460mm	313mm	1670mm	1430mm

Specifications

Defrost Mode	Refrigerant	Sound Level (dB(A))	Cooling Performance (°C)	Temperature Range (°C)	Power Supply (V/Hz)	Power (W)
Manual	HC	34	-25	-10~-25	220~240/50	270
Electrical Current (A)	Power Consumption (kW·h/24h)	Capacity (L/Cu.Ft)	Net/Gross Weight (approx)(kg)	Internal Dimension (W*D*H)(mm)		
1.4	0.95	262/9.3	88/93	480*465*1430		
Exterior Dimension (W*D*H)(mm)	Container Load (20'/40'/40'H)	Alarms				Certification
700*705*1665	21/45/45	High/Low Temperature, Remote Alarm, Power Failure, Sensor Error, High Ambient Temperature, Door Ajar				CE