


## Haier Laboratory Refrigerator

Wide Range Temperature Adjustment  
Precise Temperature Stabilization Control  
Designed specifically for labs

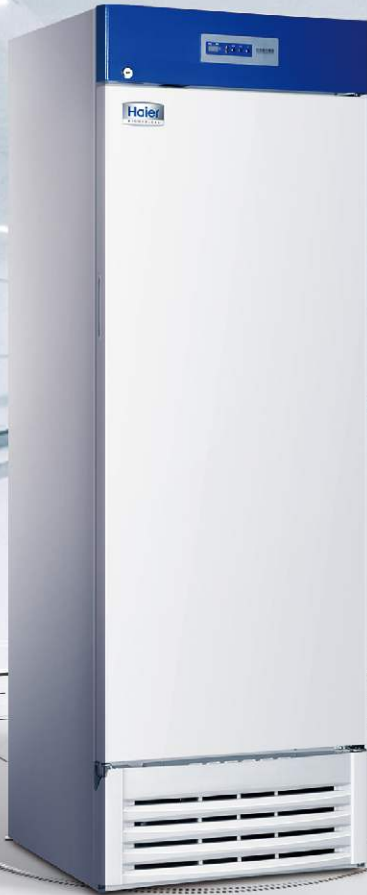

3°C

16°C

\* Temperature control is precise and stable



HLR-198F/310F

- **Wide Range Temperature Adjustment:**  
3-16°C wide temperature range designed to adjust with one key for better storing suitability.
- **Precise Temperature Stabilization Control:**  
Internal temperature variation is within 3°C to professionally keep temperature uniformity.
- **Energy Saving and Environmentally Friendly:**  
Hydro carbon cooling system and energy saving fan motor provide environment friendliness and energy saving.
- **Multiple Alarms:**  
Include high and low temperature alarms, door open alarm, etc. The refrigerator can be remotely monitored and will alert users rapidly.



### Haier Medical and Laboratory Products Co.,Ltd

Address: Brand Building, Haier Industrial Park, No.1  
Haier Road, Qingdao, 266101, P.R.China  
Tel: +86-0532-88936011 / 5955  
Website: www.haiermedical.com

**Main Uses** ▶ To store chemical and experiment reagents and consumables for institutional laboratories, research facilities, etc.

### Product Advantages



**3-16°C wide temperature range with precise control to prevent freezing**  
3-16°C wide range temperature adjustment. Intelligent temperature control system keeps internal temperature variation within 3°C for great temperature uniformity.



**One Key Search for Convenient Checking**  
Historical temperature data can be saved and checked by pressing only one key for convenient laboratory inspection.



**Dual Lock Design with Multiple Alarm System**  
The door comes with a mortise lock. A padlock can also be installed. The system has high and low temperature, open door, and sensor error alerting functions for better safety.

### Ergonomic Design



**HC Cooling Saves Half Energy**  
The hydro carbon cooling system has twice the cooling efficiency and saves half the energy.



**Low Noise**  
Ultra quiet fan motor with noise cancelling technology makes the entire refrigerator's noise lower than 40 decibels for a quiet and comfortable lab environment.



**Glass Shelves**  
Tempered glass shelves can store over 40kg of various sized reagent bottles.



**Automatic Door Closing Function**  
Provides small angle automatic door closing function.



Cabinet Structure



Glass Shelve



Fan/Interior Light



Display Board

### Product Specifications

Model	Voltage	Power	Internal Temperature Range	Exterior Dimensions (W*D*H)	Interior Dimensions (W*D*H)	Capacity	Shelves	Net Weight
HLR-198F	220V/50Hz	195W	3°C-16°C	530*570*1565	455*445*900	198L	4	54kg
HLR-310F	220V/50Hz	270W	3°C-16°C	605*600*1870	560*500*1290	310L	4	77kg